


# Spinsolve Phosphorus

## Phosphorus-31 benchtop NMR

Phosphorus is commonly found in many organic compounds, for example in biological membranes or DNA. The  $^{31}\text{P}$  nucleus has a 100% natural isotopic abundance and a large chemical shift range, making it one of the most commonly used nuclei in biological NMR.

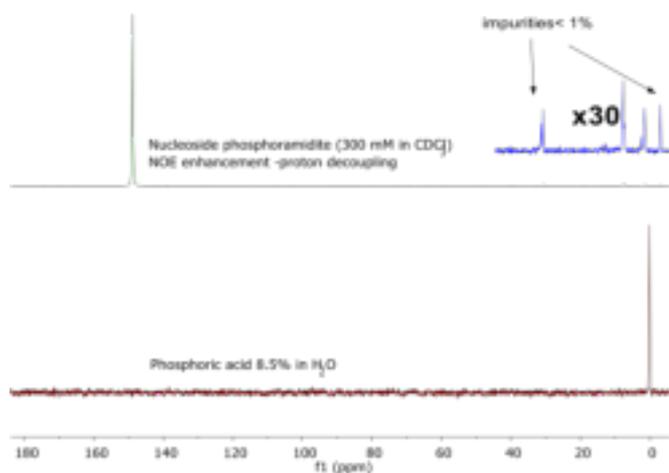

### Who


- Medicinal Chemists
- Pharmaceutical Chemists
- Structure Elucidation
- Purity Assessment

### Why


- No cryogenics
- Fast
- Convenient
- Low cost
- Accessible
- Robust
- Low maintenance
- Easy to operate
- Exceptional performance

## $^{31}\text{P}$ NMR Spectroscopy



Example NMR Spectra of Phosphoramidites using Spinsolve Phosphorus.

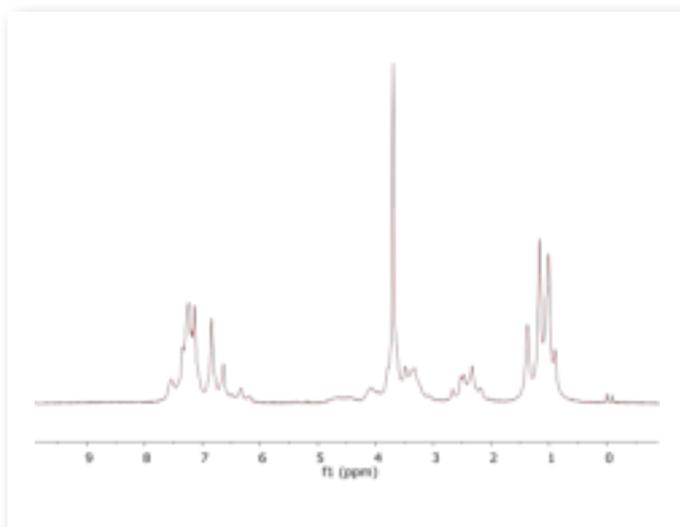

dT Phosphoramidite

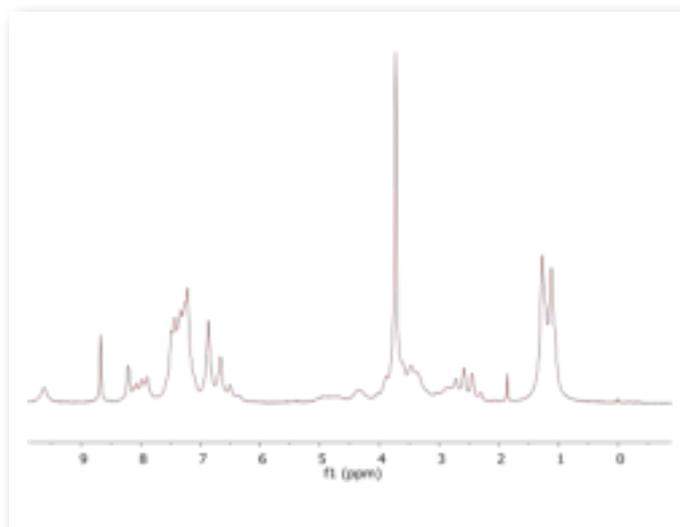

dA(bz) Phosphoramidite

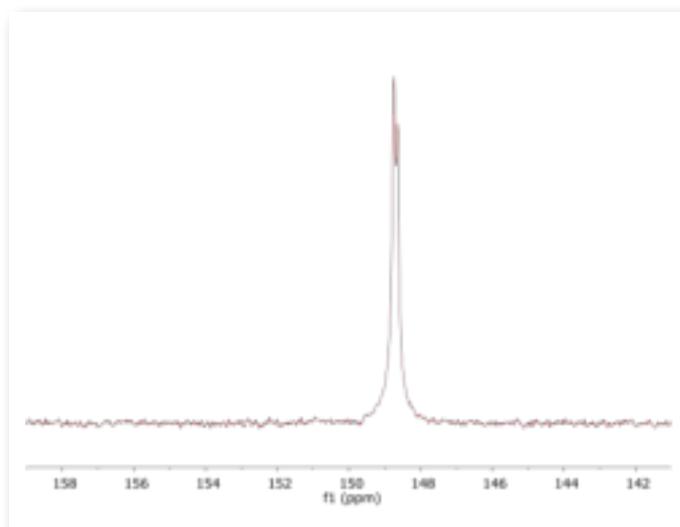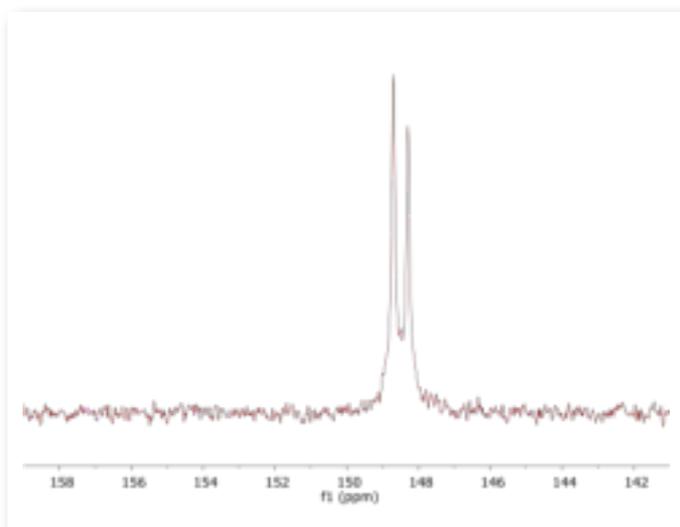


$^1\text{H}$  (top) and  $^{31}\text{P}$  (bottom) spectra of two phosphoramidite solutions (200 mM in  $\text{CDCl}_3$ ). These molecules play an important role in DNA synthesis. The measurement time for the  $^{31}\text{P}$  spectra was 7 minutes.

# Software

## The power of simplicity

- Uncomplicated one-button operation
- Simple intuitive graphical interface
- Automated and easy to use
- Minimal user controlled parameters
- Traditional NMR complexities are hidden

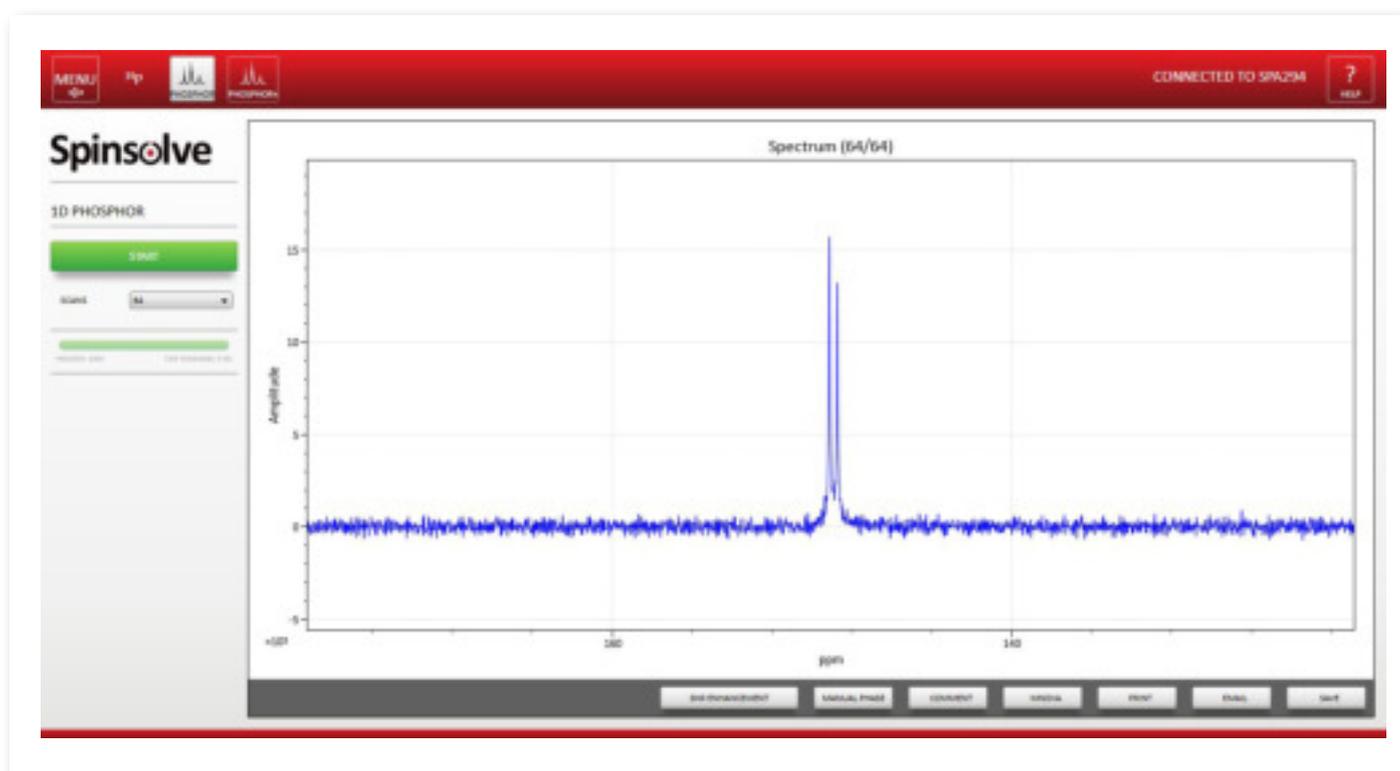

## Features

- 1D  $^1\text{H}$ ,  $^{19}\text{F}$  and  $^{31}\text{P}$  experiments
- Standard 5 mm NMR tubes
- 2D COSY and homonuclear j-resolved spectroscopy
- $T_1$  and  $T_2$  relaxation experiments
- Composite pulse decoupling

# Spinsolve Phosphorus



## Specifications

- Frequency: 42.5 MHz Proton, 17.4 MHz Phosphorus
- $^1\text{H}$  Resolution: 50% linewidth  $< 0.5$  Hz
- $^1\text{H}$  Lineshape: 0.55% linewidth  $< 20$  Hz
- Dimensions: 58 x 43 x 40 cm
- Weight: 55 kg
- Magnet: Permanent and cryogen free
- Stray field:  $< 2$  G all around system

## Other Spinsolve products

### Spinsolve<sup>®</sup>



- $^1\text{H}$  and  $^{19}\text{F}$  nuclei
- Relaxation time experiments
- 2D COSY and JRES
- Reaction monitoring

### Spinsolve<sup>®</sup> Carbon



- $^1\text{H}$ ,  $^{19}\text{F}$  and  $^{13}\text{C}$  nuclei
- Spectral editing with DEPT
- Composite pulse decoupling
- 2D heteronuclear correlation experiments HECTOR, HMQC, NMBC, HSQC

Contact us now for a quote, to request a demo or to measure your samples

Email: [sales@magritek.com](mailto:sales@magritek.com)

Website: [www.magritek.com/contact-us](http://www.magritek.com/contact-us)

GERMANY +49 241 9278 7270

UNITED KINGDOM +44 7468 529 615

UNITED STATES +1 855 667 6835

NEW ZEALAND +64 4 477 7096

For a complete listing of our global offices and distributors visit: [www.magritek.com/about-us/distributors/](http://www.magritek.com/about-us/distributors/)